

MixStar® Cartridges & Accessories

- 999537 Honigum® Heavy (1-380ml Cartridge plus 10 MixStar Tips)
- 999568 Honigum® Mono (1-380ml Cartridge plus 10 MixStar Tips)
- 999638 Honigum® QuadFast® Heavy (1-380ml Cartridge plus 10 MixStar Tips)
- 999374 Honigum® MixStar Putty (1-380ml Cartridge plus 10 MixStar Tips)
- 999693 StatusBlue® (1-380ml Cartridge plus 10 MixStar Tips)
- 999482 MixStar® eMotion Machine
- 999310 MixStar® Foot Pedal (Optional)
- 999567 MixStar® Mixing Tips (Bag of 50)

50ml Automix Cartridges & Accessories

- 999766 Honigum® Heavy (4-50ml Cartridges, 12 Automix (Green) and Intra-Oral Tips)
- 999830 Honigum® Mono (4-50ml Cartridges, 12 Automix (Green) and Intra-Oral Tips)
- 999838 Honigum® QuadFast® Heavy (4-50ml Cartridges, 12 Automix (Green) and Intra-Oral Tips)
- 999826 StatusBlue® Economy Tub (7 +1 FREE 50ml Cartridges, 16 Automix Tips)
- 110411 StatusBlue® Type 50 Dispensing Gun (10:1)
- 999400 Honigum® 50ml Tips
- 999797 StatusBlue® 50ml Tips
- 999490 Honigum® Mono Intra-Oral Tips

Tray Adhesive

- 999394 Tray Adhesive

25ml Automix Cartridges & Accessories

- 999539 Honigum® Light Body (4-25ml Cartridges, 16 Automix and Intra-Oral Tips)
- 999639 Honigum® QuadFast® Light (4-25ml Cartridges, 16 Automix and Intra-Oral Tips)
- 110253 Type 25 Dispensing Gun
- 110424 Automix and Intra-Oral Tips (35 of Each)

O-Bite™ (VPS)

- 999765 O-Bite (4-50ml Cartridges and 15 Tips)
- 999507 Automix Dispensing Gun (Type 50 1:1)
- 999202 Automix Tips (50, Pink 1:1)

LuxaBite® Bisacryl

- 110560 LuxaBite® Kit (1-50ml Cartridge, 15 Automix Tips)
- 110609 LuxaBite® Automix Tips (25 LuxaBite Tips)
- 110411 LuxaBite® Type 50 Dispensing Gun (10:1)

Honigum® Hand Mix Putty

- 999366 Rigid Fast (1-300ml Rigid Fast Base, 1-300 ml Rigid Fast Catalyst, 2 Accessory Spoons)
- 999365 Soft Fast (1-300ml Soft Fast Base, 1-300 ml Soft Fast Catalyst, 2 Accessory Spoons)

Honigum® Impression System

■ The complete impression material system.



Advantages

- A comprehensive system for all your impression needs
- Quality dental materials that provide exceptional ease and accuracy
- A choice of advanced Automix dispensing options



For more information, call us at **1-800-662-6383**. To find out more about all of our innovative dental products, visit our website at www.dmg-america.com

DMG America

242 South Dean Street
Englewood, New Jersey 07631
P: 800-662-6383 • 201-894-5505
F: 201-894-0213 • www.dmg-america.com



▲ Honigum ▲ LuxaBite ▲ StatusBlue ▲ MixStar ▲ O-Bite

The Honigum® Impression System:

For Amazingly Accurate and Easy Impressions

From high-detail Honigum tray and wash impression materials, to innovative StatusBlue Automix A-Silicone alginate replacement and superior O-Bite and LuxaBite registration materials, DMG America has everything you need to create fast, easy, and the most accurate impressions. Our superior materials and dispensing options make it easy to save time and get great results every time, eliminating costly mistakes and re-makes.

Honigum® Impression Materials: Capture Incredible Detail Every Time

The secret to Honigum's superior detail lies in its unique chemistry. Honigum is a hydrophilic poly vinylsiloxane impression material with a difference – a patented microcrystalline wax matrix that lowers its viscosity under bite pressure to flow into the most difficult to reach subgingival areas. This allows Honigum to capture details other impression materials can't! Honigum also has a "snap set" that prevents distortion, yet is easily removable. Honigum's pleasant honey aroma and neutral taste also help to make the material extremely patient-friendly.


Honigum® Tray Materials

As a tray material, Honigum is perfect for crown and bridge, inlay/onlay, implant and full mouth restorations. Honigum comes in 3 patented non-slumping formulations that are available in either 380ml cartridges for the MixStar eMotion machine (that will also fit the Pentamix® Machine), or 50 ml Automix cartridges for use with Type 50 1:1 dispensing guns. Honigum tray materials are used in conjunction with Honigum wash materials.

Pentamix® is a registered trademark of 3M ESPE

Ordinary Materials' Flow

Conventional impression materials do not lower their viscosity under bite pressure, so they do not exhibit the flow necessary to capture minute subgingival details.

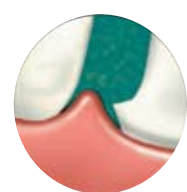

Honigum's Direct Flow
Honigum lowers its viscosity under bite pressure, reaching into even the finest subgingival crevices to capture amazing detail.

Honigum® Wash Materials

Honigum wash materials provide increased impression detail when used in conjunction with Honigum tray materials. All offer Honigum's unique patented formulation and are available in convenient 25ml Automix cartridges for the Type 25 dispensing guns.



Easy direct placement of Honigum Light.

DISPENSING KEY:

- MixStar Machine
- Type 25ml Automix
- 1:1 Type 50ml Automix
- 10:1 Type 50ml Automix
- Handmix



The contrasting colors of Honigum Heavy and Honigum Light capture and show details better.

Honigum® Heavy Body

Description A bright blue tray material
Techniques Perfect for double-phase impressions
Setting Time 3.15 min.



Honigum® Mono

Description Universal, medium viscosity, green tray and wash material
Techniques For an exact single-phase technique
Setting Time 2.30 Automix
3.15 MixStar



Honigum® QuadFast® Heavy

Description Fast-setting, teal tray material
Techniques For fast double-phase quadrant and single-unit impressions
Setting Time 2.0 min.



Honigum® Light Body

Description High-contrast bright yellow wash material
Techniques For one-step and two-step, plus double-phase quadrant and single-unit impressions
Setting Time 3.30 min.



Honigum® QuadFast® Light

Description Bright yellow-orange wash material
Techniques For one-step and two-step, plus double-phase quadrant and single-unit impressions
Setting Time 2.0 min.



The MixStar® eMotion Machine

MixStar eMotion is an efficient, dynamic mixing machine that provides clean, easy and hands-free mixing and dispensing of pre-filled 380ml cartridges. MixStar eMotion includes two built-in times indicating working and setting times for up to six impression materials.



Bite Registration Materials

Now you have a choice of superior rigid bite registrations materials: O-Bite and bisacryl LuxaBite. All are ideal for capturing all types of bite relationships, easy to trim and are available in various dispensing methods. O-Bite is easily trimmed with a scalpel and LuxaBite is easily trimmed with a high-speed burr.



Exceptional occlusal accuracy is achieved with our proprietary vinyl polysiloxane matrix.

StatusBlue® Situation Impression Material

StatusBlue is an innovative A-Silicone impression material that's specifically designed as a high efficiency alternative to traditional alginates. Use it for temporaries, model cast dentures, orthodontic models and opposing jaw models. StatusBlue can be cold disinfected and re-poured as needed for up to 2 weeks. Available in 380ml cartridges for the MixStar eMotion machine (also fits the Pentamix® machine) or 50ml Automix cartridges for use with Type 50 10:1 dispensing guns. Our convenient Automix options give you easy mixing and dispensing with no hand mixing, no mess, and no clean-up.

MixStar® Putty

Description Dynamically mixed true putty
Techniques The first VPS impression material with a putty consistency available in a MixStar cartridge
Setting Time 2.0 min.



Honigum® Putty Soft Fast

Description VPS impression material, putty consistency
Techniques Utilized in the double-mix technique in combination with Honigum Light-Body
Setting Time 2.0 min.



Honigum® Putty Rigid Fast

Description VPS impression material, putty consistency
Techniques Designed for the correction technique
Setting Time 2.0 min.



O-Bite™ (VPS)

Description A rigid-setting orange colored and scented material
Techniques Captures accurate and stable bite registrations, easily trimmed with scalpel
Setting Time 1.30 min.



LuxaBite® Bisacryl

Description An ultra rigid-setting blue composite bite registration material
Techniques Captures accurate and stable bite registrations
Setting Time 2.00-2.30 min. after start of mix
4 min. after removal from mouth



StatusBlue® Situation Impression Material

Description A soft flowable blue tray material
Techniques Excellent for all types of situation impressions
Setting Time 1.45 minutes

