


ZEKRYA GINGIVAL PROTECTOR



SWISS MADE



The Zekrya Gingival Protector provides protection for your patient and convenience for you. With the Zekrya Gingival Protector, you can retract soft tissue to protect gingiva from rotary instruments during cutting and finishing. Zekrya features an exclusive crescent-shaped tip on an adjustable ball joint that is attached to a metal handle. The tip can be rotated to an exact angle that matches the tooth's facial surface to achieve gingival fit. The protector is ideal for veneer preparations, finishing porcelain or resin veneer margins, class V (facial) subgingival caries removal, checking marginal fit of crowns, and checking loose bridges.

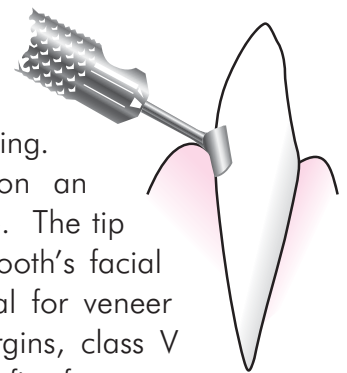

The Zekrya Gingival Protector prevents injury to the gingival tissue and causes less discomfort to the patient. Insertion of a restoration is possible without waiting for healing. When the Zekrya protector is used, there is no blood to interfere with cementation of the veneer or the polymerization of a composite. Zekrya increases direct visibility by offering more accurate and controlled preparation and interpretation of margins.

Metal protector tips are autoclavable to prevent cross-contamination. Tips come in two sizes to fit all teeth.

Product Number	Description
61096	Zekrya Kit: (10 Small Arcs, 10 Large Arcs & 1 Gingival Protector)
61097	Bag of 10 Small Arcs
61098	Bag of 10 Large Arcs

 **ZenithDental™**
Enhancing dentistry with education and technology®

242 South Dean St., Englewood, NJ 07631
(201)-894-5500 • 1-800-662-6383
Fax: 201-816-1812 • www.zenithdental.com

ZEKRYA GINGIVAL PROTECTOR



SWISS MADE

THE ZEKRYA GINGIVAL PROTECTOR ALLOWS:

- Retraction of the gum without injury-risk
- Complete protection of the gum during tooth preparation
- Good visibility

TECHNICAL FEATURES:

- Gingival protectors are sterilizable and replaceable
- Swivel joint allows an easy adjustment to any desired position
- The flexible arc of the protector adapts to fit all tooth necks
- Retracts soft tissue while protecting gingiva during rotary instrumentation.

TIPS FOR USE:

- Place the protector on the mandrel in the desired position and fix it firmly
- Protect the patient's lips by placing a cotton roll between gum and lips
- Press the protector firmly against the tooth neck so that it fits
- Spray with water in order to prevent overheating while using a rotary instrument