



Luxatemp

The ultimate esthetic provisional system



Perfect Provisionals

Developed by DMG more than 20 years ago, Luxatemp was the first bisacryl provisional material introduced as an automix material in the U.S. and Canada. Today, it remains the system of choice for clinicians to produce esthetically superior, lifelike provisionals.

Whether it's award-winning Luxatemp Fluorescence's superior esthetics, Luxatemp Ultra's flexural strength and break resistance, or the proven reliability of Luxatemp Plus, the complete DMG family of provisional materials continues to deliver unmatched clinical results and premium quality to meet the needs of dentists and patients alike.

Luxatemp® Ultra

Rated #1 in flexural strength*

Luxatemp Ultra continues to deliver the unmatched esthetics synonymous with the Luxatemp name. Along with superior color stability, it includes fluorescence to create temporaries that look more like natural teeth in all lighting conditions.

FEATURES & BENEFITS

- Industry leader in flexural strength
- Unsurpassed break resistance
- Faster curing, easier cutting and handling for greater clinician control
- Fluorescence added for outstanding esthetics
- Superior stability and long-term durability

Luxatemp® Fluorescence

Naturally esthetic

- Natural fluorescence
- Lifelike provisionals look great in all lighting conditions
- High fracture and abrasion resistance
- Durable, long-lasting temporaries
- Excellent marginal adaption
- Easy to remove
- Minimal shrinkage during polymerization
- Outstanding color stability




Before Luxatemp provisionals



After Luxatemp Provisionals

Photos courtesy of:
Dr. Jon Marashi, DDS, San Clemente, CA

FLEXURAL STRENGTH

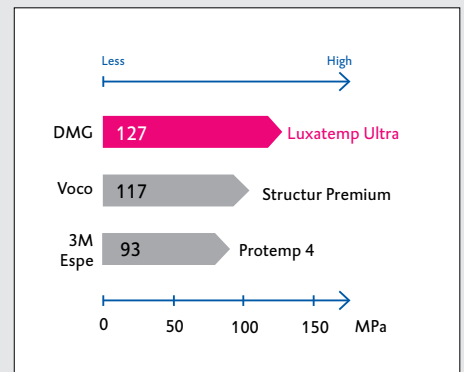

Measurements carried out at Erlangen University Clinic showed outstanding flexural strength values for Luxatemp Ultra.¹

LUXATEMP® ULTRA PHYSICAL PROPERTIES

Chemistry: bisacryl

Mixing ratio (base/catalyst): 10:1 (vol.)

Compressive strength: 376 MPa

Transverse strength: 100 MPa

Diametrical tensile strength: 40 MPa

Barcol hardness: 40

Water sorption: 20µ/mm³

Temp. rise during polymerization: 39°C/102°F

* "Flexural Strength of Provisional Restoratives," a comparative study; Perry R., et al

Everybody Has a Favorite

Luxatemp® Automix Plus

- Durable and long-lasting provisionals
- Exceptional precision fit
- High biocompatibility
- High break resistance
- Excellent flexural strength
- Low shrinkage during polymerization

Luxatemp® Solar

- “Cure-on-demand” for extended working time
- Easy removal during elastic phase
- High final hardness
- Low shrinkage during polymerization
- Excellent polishing properties

Complete Your Provisional System



TempoCem® ID Temporary Cement

Invisible under restorations, yet detectable for excess removal

- Beautiful esthetics thanks to optimal transparency
- Easy one-piece removal of excess
- Optimal adhesion with easy removability
- Eugenol free
- Excellent ease of use, no sticking to instruments
- Flows and mixes easily
- Perfect for long-term temporaries (up to 6 mos.)



LuxaFlow™ Ultra/LuxaFlow™

Ideal for lifelike repairs

- Optimum flow characteristics
- Highly polishable
- Non-slumping formula
- Light cured
- Radiopaque

LuxaGlaze®

For a high-gloss finish

- High gloss shine without polishing
- Greater protection, stain-resistant
- Exhibits life-like appearance
- Enhances temporary esthetics
- Quick and easy application
- Light cured



Luxatemp Ultra material being dispensed into a StatusBlue®* matrix.



Removing the temporary from the matrix.



Refining the temporary with a bur.



A highly esthetic, natural-looking temporary end result.

*StatusBlue alginate replacement

Photos courtesy of:
Dr. Tony Soileau, DDS, Lafayette, LA



Ordering Information



Luxatemp Automix Plus

Item #	Product Description	Contents
110400	Luxatemp Plus Automix Intro Kit - A2	1-76g Cartridge, 15 Automix Tips and 1 Type 50 Dispensing Gun
110401	Luxatemp Plus Automix Refill—A1	1-76g Cartridge, 15 Automix Tips
110402	Luxatemp Plus Automix Refill—A2	1-76g Cartridge, 15 Automix Tips
110403	Luxatemp Plus Automix Refill—A3.5	1-76g Cartridge, 15 Automix Tips
110404	Luxatemp Plus Automix Refill—B1	1-76g Cartridge, 15 Automix Tips
110368	Luxatemp Plus Automix Refill—Bleach Light	1-76g Cartridge, 15 Automix Tips
110406	Luxatemp Plus Automix Bulk Pack—A2	5-76g Cartridges and 75 Automix Tips

Luxatemp Fluorescence

110585	Luxatemp Fluorescence Automix Refill—A1	1-76g Cartridge, 15 Automix Tips
110586	Luxatemp Fluorescence Automix Refill—A2	1-76g Cartridge, 15 Automix Tips
110710	Luxatemp Fluorescence Automix Refill—A3	1-76g Cartridge, 15 Automix Tips
110587	Luxatemp Fluorescence Automix Refill—A3.5	1-76g Cartridge, 15 Automix Tips
110588	Luxatemp Fluorescence Automix Refill—B1	1-76g Cartridge, 15 Automix Tips
110589	Luxatemp Fluorescence Automix Refill—BL	1-76g Cartridge, 15 Automix Tips
110719	Luxatemp Fluorescence Smartmix—A2	1-15g Syringe, 10 Smartmix Tips
110722	Luxatemp Fluorescence Smartmix—B1	1-15g Syringe, 10 Smartmix Tips
110723	Luxatemp Fluorescence Smartmix—BL	1-15g Syringe, 10 Smartmix Tips

Luxatemp Ultra

110906	Luxatemp Ultra Automix Refill—A1	1-76g Cartridge, 15 Automix Tips
110907	Luxatemp Ultra Automix Refill—A2	1-76g Cartridge, 15 Automix Tips
110908	Luxatemp Ultra Automix Refill—A3	1-76g Cartridge, 15 Automix Tips
110909	Luxatemp Ultra Automix Refill—A3.5	1-76g Cartridge, 15 Automix Tips
110910	Luxatemp Ultra Automix Refill—B1	1-76g Cartridge, 15 Automix Tips
110911	Luxatemp Ultra Automix Refill—BL	1-76g Cartridge, 15 Automix Tips
110901	Luxatemp Ultra Smartmix Refill—A2	1-15g Syringe, 10 Automix Tips
110904	Luxatemp Ultra Smartmix Refill—B1	1-15g Syringe, 10 Automix Tips
110905	Luxatemp Ultra Smartmix Refill—BL	1-15g Syringe, 10 Automix Tips

Luxatemp Automix Solar

110360	Luxatemp Automix Solar Refill—A2	1-76g Cartridge, 15 Automix Tips
--------	----------------------------------	----------------------------------

LuxaFlow/Luxaflow Ultra

211751	LuxaFlow Intro Kit	1-1.5g Syringe each of A1, A2, A3.5 and B1, plus 20 Syringe Tips
211753	LuxaFlow A1 Refill	2-1.5g Syringe and 10 Syringe Tips
211754	LuxaFlow A2 Refill	2-1.5g Syringe and 10 Syringe Tips
211755	LuxaFlow A3.5 Refill	2-1.5g Syringe and 10 Syringe Tips
211756	LuxaFlow B1 Refill	2-1.5g Syringe and 10 Syringe Tips
224000	LuxaFlow Ultra Intro Kit	1-1.5g Syringe each of A1, A2, A3.5 and B1, plus 20 Syringe Tips
224001	LuxaFlow Ultra A1 Refill	2-1.5g Syringe and 10 Syringe Tips
224002	LuxaFlow Ultra A2 Refill	2-1.5g Syringe and 10 Syringe Tips
224003	LuxaFlow Ultra A3 Refill	2-1.5g Syringe and 10 Syringe Tips
224004	LuxaFlow Ultra A3.5 Refill	2-1.5g Syringe and 10 Syringe Tips
224005	LuxaFlow Ultra B1 Refill	2-1.5g Syringe and 10 Syringe Tips

LuxaGlaze

212075	LuxaGlaze Kit	5ml Dropper Bottle, Mixing Pallet and 25 Applicator Brushes
--------	---------------	---

TempoCem®ID, TempoCem, TempoCem NE

212075	TempoCemID SmartMix	1-5ml Syringe and 10 Smartmix Tips
110300	TempoCem Automix Refill Kit	1-25ml Cartridge and 40 Automix Tips
110302	TempoCemNE Automix Refill Kit	1-25ml Cartridge and 40 Automix Tips
212104	TempoCemNE Smartmix	2-11g Syringes and 10 Smartmix Tips
212129	TempoCem Smartmix	2-11g Syringes and 10 Smartmix Tips

Accessories

110411	Type 50 10:1/4:1 Applicator Gun	1 Gun
110253	Type 25 1:1 Applicator Gun	1 Gun
110351	Luxatemp Automix Mixing Tips	Bag of 25
110991	Luxatemp Smartmix Orange-Green Tips	Bag of 25
211759	LuxaFlow Ultra/LuxaFlow Syringe Tips	Bag of 20
212040	TempoCem Smartmix Mixing Tips	Bag of 50
110266	TempoCem Automix Mixing Tips	Bag of 50
110728	TempoCemID Smartmix Tips	Bag of 25

“I’m a longtime Luxatemp user, but Luxatemp Ultra takes your restorations to the next level with increased polishability and overall esthetics.”

- Dr. Courtney Lavigne, DMD, FAGD

DMG America, LLC
 65 Challenger Road, Suite 340, Ridgefield Park, NJ 07660
 tel: 800.662.6383 email: info@dmg-america.com
 web: www.dmg-america.com

