

00:01:00

ONE MINUTE PRODUCT MEETING


Ionosit Baseliner

Ionosit Baseliner is a single component, micro-expansion agent that prevents micro-gaps, the leading cause of microleakage, before they happen.

Capable of controlled expansion thanks to its innovative chemistry, Ionosit Baseliner also forms a connection to dentin. This creates a so-called reaction zone. The crucial advantage: possible gaps created by polymerization shrinkage of the composite material are negated and stress on the dentin interface is reduced.

Plus, it releases fluoride and zinc ions that that provides protection at the junction of the tooth surface and the restoration.

Ionosit Baseliner is a base material for composite and amalgam restorations that can be used to block out undercuts in inlay/onlay restorations.

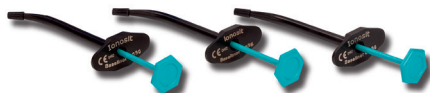

FEATURES & BENEFITS

- › Micro-expansion agent
- › Bonds to dentin and composite and protects pulp
- › Fluoride release: antibacterial properties
- › Bacteriostatic
- › Light cured: plenty of working time
- › Radiopaque
- › No complicated hand-mixing

APPLICATIONS

- › Baseliner for composite restorations
- › Baseliner for amalgam restorations

ITEM # & DESCRIPTION

210911 | Ionosit Baseliner Office Pack

223002 | Ionosit Baseliner Syringe Pack

LOOKING FOR MORE INFO?

Visit www.dmg-america.com/products.